
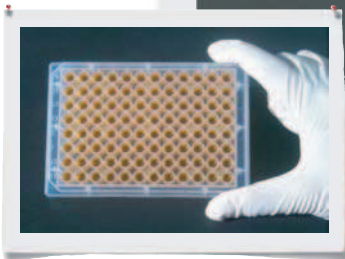


## PolyPlex®

High-throughput, oligo synthesis in a 96-well format.

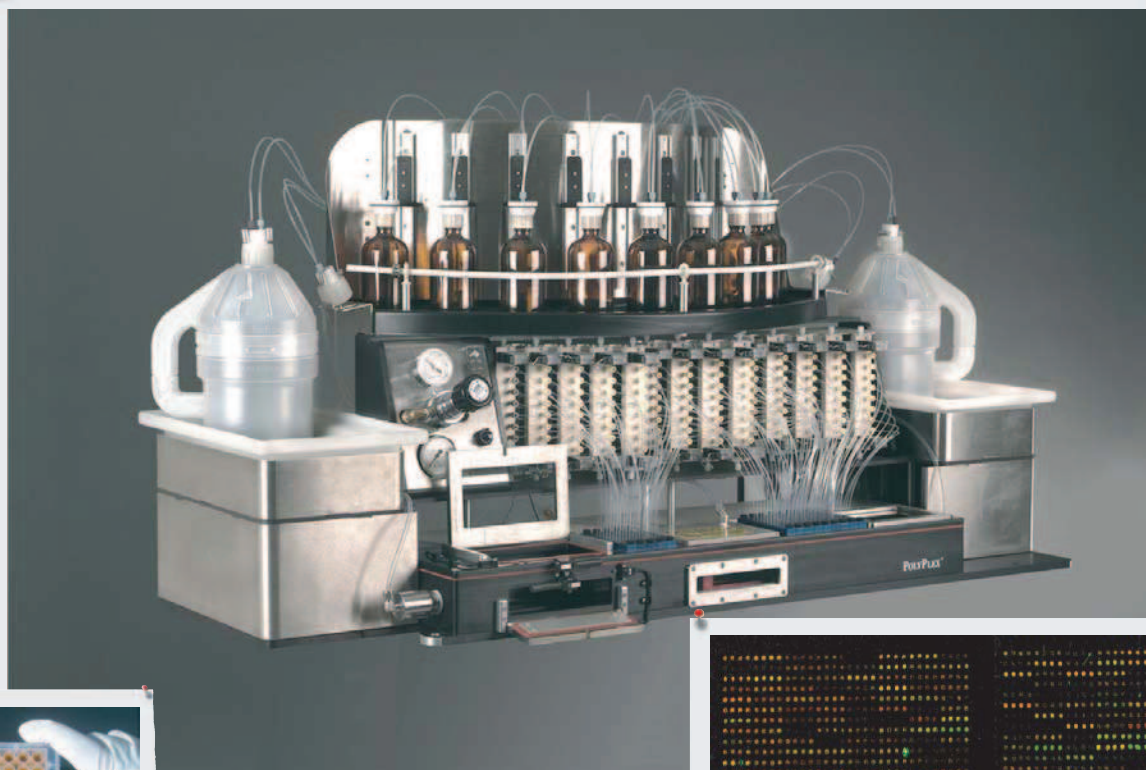

### Trityl Collection at Any Cycle

The PolyPlex is a fast, cost-effective, 96-well oligonucleotide synthesizer capable of generating a full plate of 20-mers in less than 3 hours. The second generation PolyPlex features up to eight amidite positions to meet the ever-growing demand for synthesizing specialty oligos.

### Oligo Array Hybridization

PolyPlex Oligos on  
OmniGrid Aldehyde Slides

Data courtesy of Roger Kroes at  
Dept. of Biomedical Engineering,  
Falk Center for Molecular Therapeutics,  
Northwestern University

Highly reliable valves, with low dead volumes and accurate dispensing, deliver low cost synthesis at **less than \$0.10 per base**. Different scales are available, up to 1  $\mu\text{mol}$ , using cartridges. Trityl collection is available after any base addition providing in situ QC. Each PolyPlex instrument

includes one year of technical and onsite-service support. Validated protocols are provided along with **powerful PolyPlex software**. Flexible importing of oligo sequences and synthesis in 96-well format are ideal features for seamless integration with the product line of **OmniGrid** microarrays, slides and reagents.



### Key Benefits

- Up to 4 auxiliary positions for special chemistries
- Less than \$0.10 per base
- Convenient 96-well format
- User-friendly software with flexible sequence loading

DNA Shearing

Colony Picking

Cell Growth

DNA/RNA Synthesis

Microarraying

Hybridization


Scanning

# PolyPlex<sup>®</sup>


## Design Specifications

<b>Dimensions</b>	60" W x 21" D x 31" H (1.52m x 0.53m x 0.79m)
<b>Cost per Synthesis</b>	\$0.086 per base for a 20-mer synthesis; \$0.078 per base for a 70-mer synthesis. Includes all costs needed to make ready-to-use crude oligos at 50 nmol scale.
<b>Reagent Positions</b>	Up to 15 positions available to accommodate for specialty chemistries
<b>Capacity per Run</b>	96 wells; up to 1 µmol scale each well.
<b>Argon Consumption</b>	55 std. ft <sup>3</sup> for a full plate of 20-mers, or 6 plates per standard cylinder (type 300 cylinder).


## Data and Software

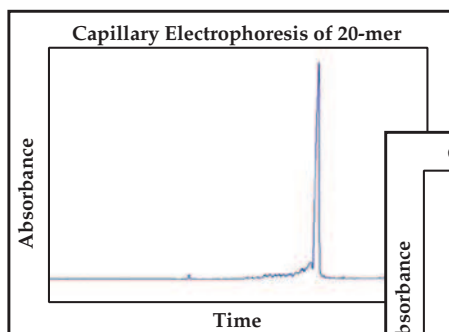


20-mer prepared with the PolyPlex in 2.7 hours.  
Coupling Efficiency- 98.8%  
CGATGCTATGGATGGATG

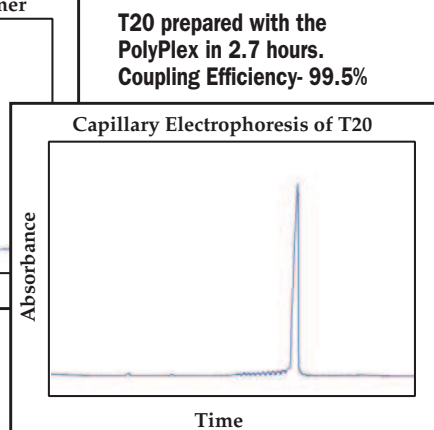

T20 prepared with the PolyPlex in 2.7 hours.  
Coupling Efficiency- 99.5%

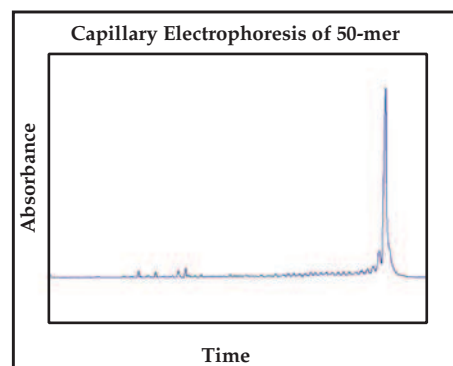


50-mer prepared with the PolyPlex in 7 hours.  
Coupling Efficiency- 99.4%  
GAATTCGAGACCACTGGTGCCTACCCCGCATCCACTCCCTCTGGC

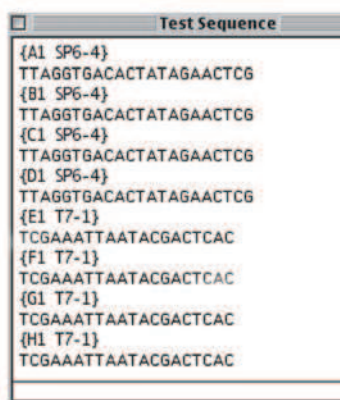

Sequence input text file format easily integrates into any database environment.

PolyPlex software provides thorough synthesis monitoring and flexible hardware control.

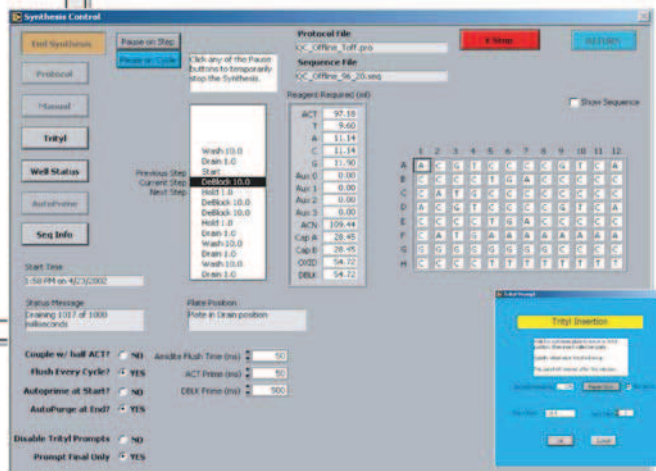

## Worldwide Headquarters

Digilab, Inc.  
84 October Hill Road  
Holliston, MA 01746  
USA

Phone: (508) 893-3130  
Toll Free: (800) 935-8007  
Fax: (508) 893-8011  
E-Mail: info@digilabglobal.com

# DIGILAB<sup>®</sup>